
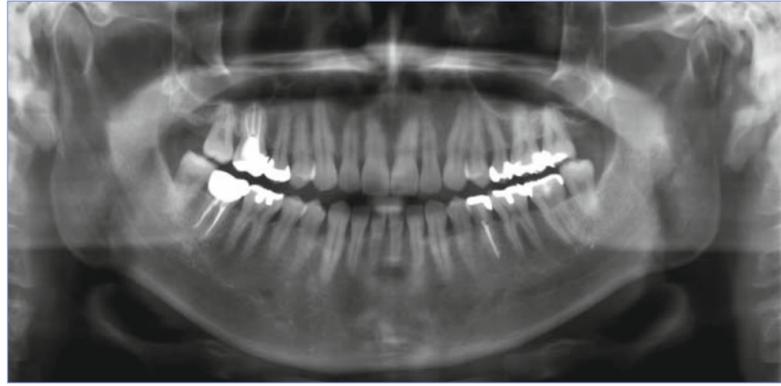


## High-quality images


### Standard panoramic

Clear, sharp images with a wide image layer. The thick, specifically designed image layer accommodates all possible variations of dental arch shapes and sizes to produce extremely clear and sharp images.

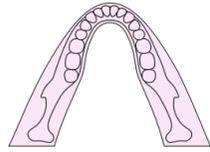



### Pedodontic panoramic at reduced radiation

For examinations of children or people with small jaws, the rotation range of the arm is diminished, thereby further reducing radiation exposure.

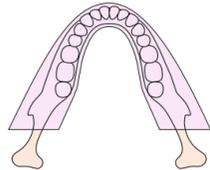

### TMJ quadruple exposures: For open and closed condyle views

Four separate shots of condyle fit onto one image. In each case, two sets of rotation for left and right condyle head for open and closed condition.

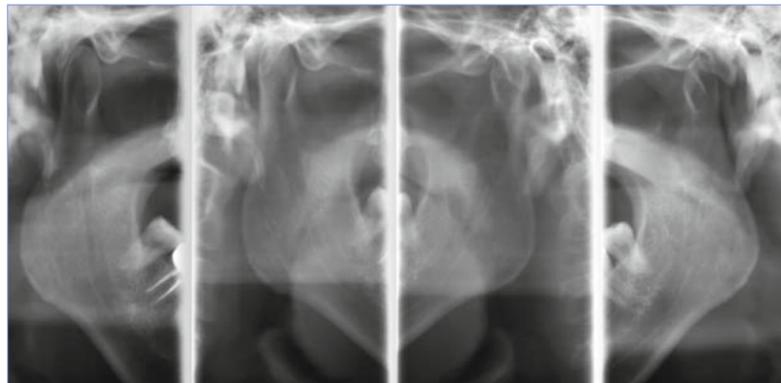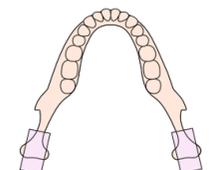

Diagnostic and Imaging Equipment

Treatment Units

Handpieces and Instruments

Endodontic Systems

Laser Equipment

Laboratory Devices



Veraview IC5 HD

High definition, digital imaging excellence



Distributed by  
**J. Morita Europe GmbH**  
 Justus-von-Liebig-Str. 27a  
 63128 Dietzenbach  
 Germany  
 T +49. 6074. 836 0, F +49. 6074. 836 299  
 www.morita.com/europe

Developed and manufactured by  
**J. Morita Mfg. Corporation**  
 680, Higashihama Minami-cho, Fushimi-ku,  
 Kyoto, 612-8533  
 Japan  
 T +81. 75. 611 2141, F +81. 75. 622 4595  
 www.morita.com

Subject to technical changes and errors. JME EN PUB 47005 0213 \*3.

Thinking ahead. Focused on life.





Thinking ahead. Focused on life.



## IC5 HD – Super High-Speed with High Definition Clarity

**For dental radiology image clarity and detail are of the essence to make a correct diagnosis. Veraview IC5 HD is the next generation panoramic X-ray machine that will give you digital exposures with outstanding image clarity, contrast, and resolution at a super high speed. The super high speed exposure takes only 5.5 seconds and therefore gives the patient a very low dose. The super high speed mode in combination with the new high definition mode will fill all your panoramic imaging needs.**

Equipped with cutting-edge technology, like automatic exposure control and automatic image enhancement, and efficient positioning, dental imaging has never been easier.

Veraview IC5 HD forms a harmonious balance between user-friendly handling and technological innovation, allowing you to conduct examinations efficiently without any sacrifice.

### Features at a glance:

#### **Industry-best exposure time\***

Extremely fast 5.5 second panoramic with significantly lower dosage.

#### **High definition**

New and improved definition of 96  $\mu\text{m}$  (10-second scan) for the ultimate high-definition image.

#### **Class-leading clarity**

DDAE (Digital Direct Auto Exposure) and AIE (Auto Image Enhancement) functions provide sharp and clear images.

#### **Quick and easy patient positioning**

With integrated triple laser beam alignment.

#### **Fully automatic simplified operation**

Simply select the desired program: Panoramic, TMJ or PEDO. In addition, no parameter setting is required.

#### **Digital technology**

No film or film development required.

#### **Low power consumption**

Lightweight and compact design.


\* as of January 2010

## High-quality images

### Digital Direct Automatic Exposure (DDAE)

DDAE controls the X-ray tube current (mA) simultaneously by detecting X-rays passing through the patient. This improves the dynamic range, and, along with Automatic Exposure (AE), results in exceptionally clear images with the best possible contrast and even density. The automatic exposure level can be adjusted to meet your individual requirements.

*There is no need to set tube voltage or current. Auto exposure (DDAE) provides the optimum exposure conditions with a maximum tube current of 7.5 mA.*

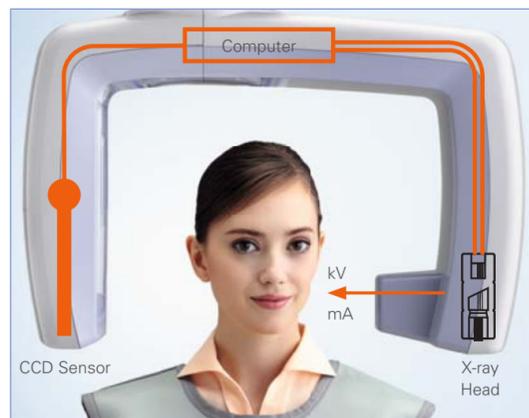

### High speed

With an industry-best exposure time\* of 5.5 seconds, the Veraview IC5 HD is not only significantly faster than conventional film-based panoramic systems, but also requires much less radiation.


High speed:  
only 5.5 seconds, 192  $\mu$ m,  
1/6\*\* X-ray exposure.

### High definition

Pixel size is reduced 25 % compared with the former model, so it produces superior images of a higher resolution.

High definition:  
10 seconds exposure, 96  $\mu$ m

\* as of January 2010  
\*\* Compared to a standard panorama image using film and made by Morita's Veraviewepocs (16 seconds, 75 kV, & 9 mA)

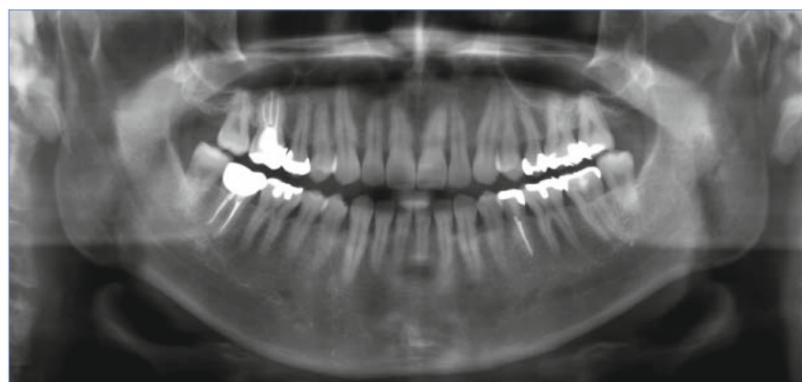

## Quick & easy patient positioning and operation

### Power-assisted lift makes patient positioning easier than ever. No need to set tube voltage and current.



#### Power-assisted movement

The motor-assisted lift system with slow-start and slow-stop function provides the easy and accurate height adjustment. The overload detector will automatically halt the elevation for the patient safety.



### Triple laser beam for easy positioning

Veraview IC5 HD utilizes three lasers that align the Frankfurt plane, the mid-sagittal plane and the image layer, making patient positioning very easy.

### Direct current for optimum efficiency and safety

Direct current produces high energy X-rays, thereby reducing the soft X-rays which are absorbed by the human body during the exposure. Also, the primary slit of Veraview IC5 HD is narrower. The image has a wider image layer, providing greater detail.



### Fully automatic X-ray loading factor

X-ray exposures require no kV or mA parameter settings. Simply select the desired program: Panorama, TMJ or PEDO.




Manually adjustable bite block/ chinrest for the patient positioning.



Chin rest cup is available for patients with TMJ disorder.

## Specifications and dimensions and patient communication features

Veraview IC5 HD · Specifications	
Trade name	Veraview IC5 HD
Model	XDP1
Input voltage	AC100 V/115 V/120 V (Ex-1) 220 V/230 V/240 V (Ex-2)
X-ray head (DC Inverter)	
Power consumption	0.85 kVA
Operating tube voltage	60-70 kV Automatic control
Operating tube current	1-7.5 mA Automatic control
Filament	Preheated
Rectification	Direct current
Cooling	Oil cooling
Tube focal area	0.5 mm
Exposure time	High speed mode: 5.5 sec./192 $\mu$ m High-definition mode: 10 sec./96 $\mu$ m
Light beams	Mid-sagittal plane, Frankfurt plane Image layer light beam
Emission switch	Deadman type
Outer dimensions	
Body dimension	W 890 x D 970 x H 2350 mm
Control box dimension	W 125 x D 60 x H 120 mm
Required installation area	0.86 m <sup>2</sup> (8.7 sq. ft.)
Weight	Approx. 110 kg (242.5 lb.)

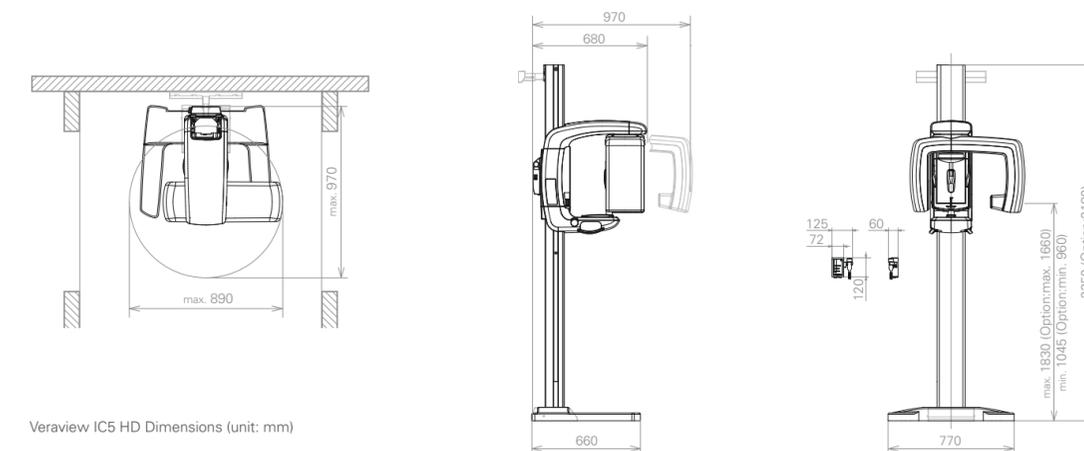

Veraview IC5 HD Dimensions (unit: mm)